



MISSOURI DEPARTMENT OF HEALTH
STATE PUBLIC HEALTH LABORATORY
DATAMASTER MAINTENANCE REPORT

BREATH ALCOHOL PROGRAM

JUN -8 2009

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.

DATAMASTER SN 201181	DATE OF INSPECTION 06/03/2009
LOCATION OF INSTRUMENT (STREET AND CITY) 120 South Ash, Nevada	TIME OF INSPECTION 1045

CHECKLIST: Place a check (✓) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> TIME AND DATE	
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34 °C ± 0.2°C) 34.0 C	
<input checked="" type="checkbox"/> CALIBRATION CHECK - Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE RECIRCULATION PUMP) <input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE <input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	

TEST 1 <input checked="" type="checkbox"/> .099	TEST 2 <input checked="" type="checkbox"/> .099	TEST 3 <input checked="" type="checkbox"/> .099
---	---	---

<input checked="" type="checkbox"/> PERFORM R.F.I. TEST (PRINTOUT ATTACHED)							
<input checked="" type="checkbox"/> NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)							
REFUSALS 2	(0-.04) 1	(.05-.09) 3	(.10-.14) 5	(.15-.19) 0	(Over .19) 0		

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary)

RepCo Marketing Inc. Lot # 08002 Exp. Date: October 13, 2010

Instrument operating satisfactorily and within specifications

INSPECTING OFFICER

SIGNATURE 	PRINT NAME Det. Jens D Barclay
TYPE II PERMIT NUMBER/EXPIRATION DATE 820113 03/26/2010	TELEPHONE NUMBER (417) 448-5152

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 08002

EXPIRATION DATE: October 13, 2010 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:

RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 08002 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed utilizing a gas chromatograph and found to contain 0.1209 gms/dl wt/vol ethyl alcohol in aqueous solution (i.e. ethanol).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of 0.100% +/-3% when heated to 34 Degrees Celsius +/-0.2 Degrees Celsius in a simulator.

The date of manufacture for this lot number is October 14, 2008. The expiration date for this lot number is October 13, 2010 at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.



Cecil B. Garner, President
RepCo Marketing, Inc.

CERTIFICATE OF ANALYSIS

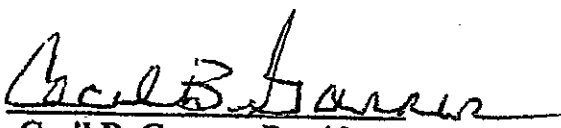
Random samples of lot number 08002 of Alcohol Certified Solution for simulator were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1209 gms/dl wt. /vol. ethyl alcohol.

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

When used in a calibrated simulator, operating at 34 c +/- .2 c, this solution will give an alcohol breath test instrument reading of .100 percent BAC +/- 2% or .002 BAC (whichever is greater).

The expiration date for this lot number is October 13, 2010 at 11:59PM.

This document is a true representation of the original Certificate of Analysis.


Cecil B. Garner, President
RepCo Marketing, Inc.

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
NEVADA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201181
06/03/09

TESTING OFFICER:
BARCLAY/JENS/D
OFFICER I.D.: 136
PERMIT NUMBER: 820113
EXPIRATION DATE: 03/26/10
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	10:47
INTERNAL STANDARD	VERIFIED	10:48
EXTERNAL STANDARD	.099	10:48
BLANK TEST	.000	10:49
EXTERNAL STANDARD	.099	10:49
BLANK TEST	.000	10:50
EXTERNAL STANDARD	.099	10:50
BLANK TEST	.000	10:51

N = 3
SIM. = .1
RMG. = .099

Operator Signature *Det James Barclay*
Printed on recycled paper with agri-based ink
CMSU 2208-02

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
NEVADA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201181
06/03/09

ARREST TIME: 10:00
SUBJECT NAME:
RFI/TEST

DOB: 06/03/09 SEX: M

STATE/D.L.: MO/1

ARRESTING OFFICER:

BARCLAY/JENS/D

OFFICER I.D.: 136

TESTING OFFICER:

BARCLAY/JENS/D

OFFICER I.D.: 136

PERMIT NUMBER: 820113

EXPIRATION DATE: 03/26/10

MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST	.000	10:53
INTERNAL STANDARD	VERIFIED	10:53
SUBJECT SAMPLE	.000	10:54
RADIO INTERFERENCE		

Operator Signature *Det James Barclay*
Printed on recycled paper with agri-based ink
CMSU 2208-02

BAC DataMaster
Evidence Ticket

STATE OF MISSOURI
NEVADA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 201181
06/03/09
10:45

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM: OKAY

HEATERS

SAMPLE CHAMBER: 49C

FLOW DETECTOR: OKAY

PUMP: OKAY

HIGH SPEED: OKAY

DETECTOR: OKAY

FILTERS: OKAY

QUARTZ STANDARD: OKAY

CALIBRATION: OKAY

PRINTER TEST

! " # \$ % & ' () * + , - . / 0 1 2 3 4 5 6 7 8 9 : ; < = > ? @ A B C D E F G
H I J K L M N O P Q R S T U V W X Y Z [\] ^ _ ` a b c d e f g h i j k l m n o
p q r s t u v w x y z { | } ~

Operator Signature *Det James Barclay*
Printed on recycled paper with agri-based ink
CMSU 2208-02

State of Missouri
DEPARTMENT OF HEALTH



P E R M I T
TYPE II



JENS DUWAIN BARCLAY

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.

Date 03/26/08

Number 820113

Expires 03/26/2010

MO 580-0771 (7-88)

Eric C. Pollock
Director of State Public Health Laboratory

[Signature]
Director, Department of Health

Lab. 4 (R7-88)