MISSOURI DEPARTMENT OF HEALTH STATE PUBLIC HEALTH LABORATORY DATAMASTER MAINTENANCE REPORT

BREATH ALCOHOL PROGRAM

JUN -7 2009

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.				
DATAMASTER SN 201302		DATE OF IN	SPECTION	
LOCATION OF INSTRUMENT (STREET AND CITY)		TIME OF IN		
501 Faraon Street, St. Joseph, Missouri	•		0806	
CHECKLIST: Place a check (1) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.				
☐ DIAGNOSTIC CHECK (PRINTOUT ATTACHED)				
☑ COMPUTER	☑ DETECTOR			
☑ PROGRAM	☑ FILTERS		·	
☑ HEATERS SAMPLE CHAMBER 49 ·c	☑ QUARTZ STANDARD	•		
☑ FLOW DETECTOR	☑ CALIBRATION		· · · · · · · · · · · · · · · · · · ·	
☑ PUMP HIGH SPEED '	☑ PRINTER			
☑ INDICATOR LIGHTS	•			
TIME AND DATE				
☑ SIMULATOR TEMPERATURE (34 -C ± 0.2-C) 34.0 •c				
☑ CALIBRATION CHECK -	•		-	
Run three tests using a standard solution. All three tests must be within $\pm 5\%$ of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE				
RECIRCULATION PUMP)		,		
☑ 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE ☐ 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE				
(ONLY ONE STANDARD IS TO BE USED PER MAINTENA				
TEST 1 .097	.098 TES	T 3 *	.098	
PERFORM R.F.I. TEST (PRINTOUT ATTACHED)				
NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)				
REFUSALS 3 (004) 5 (.0509) 3	(.1014) 5 (.1519)	8	(Over .19) 5	
List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily				
and within established limits (use other side if necessary) Instrument is operating within the limits set for by the Missouri Department of Health. Solution was				
provided by Repco Marketing. It is a .100 solution lot number 08002 expiration date of 10/13/2010				
provided by Nepco Marketing. It is a . 100 solution lot number 00002 expiration date of 10/13/2010				
· · · · · · · · · · · · · · · · · · ·		·	·	
INSPECTING OFFICER				
SIGNATURE (1)	PRINT NAME			
TYPE II PERMIT NUMBERIEXPIRATION DATE	Scott Gary Chris M	lcBane		
920051 03-06-11 820270 09-18-10	(816)271-5359	•	,	

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 08002

EXPIRATION DATE: October 13, 2010 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:

RepCo Marketing, Inc. manufactured, tested and supplied Lot Number <u>08002</u> of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed utilizing a gas chromatograph and found to contain <u>0.1209</u> gms/dl wt/vol ethyl alcohol in aqueous solution (i.e. ethanol).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of <u>0.100%</u> +/-3% when heated to 34 Degrees Celsius +/-0.2 Degrees Celsius in a simulator.

The date of manufacture for this lot number is <u>October 14, 2008</u>. The expiration date for this lot number is <u>October 13, 2010</u> at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

Cecil B. Garner, President RepCo Marking, Inc.

Form MO 01

CE THIS SIDE DOWN - THIS EDGE IN FIRST

BAC DataMaster

Evidence Ticket

ST. JOSEPH POLICE DEPT. BAC DATAMASTER SERIAL NUMBER 201302 06/06/09

RREST TIME: 07:00 JUBJECT NAME: IOB: 09/09/09 SEX: M TATE/D.L.: MO/977 RRESTING OFFICER: HJ FFICER I.D.: 865 ESTING OFFICER: **GARY/SCOTT** FFICER I.D.: 1013 'ERMIT NUMBER: 920051 XPIRATION DATE: 03/06/11 HISCELLANEOUS DATA: RFI TEST

--- BREATH ANALYSIS ---

LANK TEST .000 08:08 MTERNAL STANDARD VERIFIED 08:08 ADIO INTERFERENCE

FACE THIS SIDE DOWN - THIS EDGE IN FIRST

BAC DataMaster

Evidence Ticket

ST. JOSEPH POLICE DEPT. BAC DATAMASTER SERIAL NUMBER 201302. 06/06/09 98:96

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY PROGRAM: OKAY HEATERS SAMPLE CHAMBER: 49c FLOW DETECTOR: OKAY PUMP HIGH SPEED: OKAY DETECTOR: OKAY FILTERS: **ÖKAY** QUARTZ STANDARD: OKAY

PRINTER TEST !"##%&'<>*+,-./0123456789:;<=>?@ABCDEFG HIJKLMMOP@RSTUVWXYZ[\]^_'abcdefghijklmno pgrstuvwxyz{l}→

OKAY

I SIGNATURE

REORDER ALL SUPPLIES FROM N.P.A.S. P.O. BOX 1435; MANSFIELD, OH (44901) **OPERATOR SIGNATURE**

CALIBRATION:

Card Stock No.

WHO THE BERGING

REORDER ALL SUPPLIES FROM N.P.A.S.

THIS SIDE DOWN - THIS EDGE IN FIRST

BAC DataMaster

Evidence Ticket

ST. JOSEPH POLICE DEPT. NC DATAMASTER SERIAL NUMBER 201302 06/06/09

TING OFFICER:
#ARY/SCOTT
ICER I.D.: 1013
MIT NUMBER: 920051
IRATION DATE: 03/06/11
CELLANEOUS DATA:


--- SUPERVISOR MODE ---

NK TEST	.000	08:34 .
ERNAL STANDARD	VERIFIED	Ø8:35
ERNAL STANDARD	. 097	. 98:35
MK TEST	.000	08:36 .
ERNAL STANDARD	.098	98:3 <u>6</u>
NK TEST	.900	Ø8:36
ERNAL STANDARD	.098	Ø8:37.
NK TEST	. 999	Ø8:37

. = .1 . = .0976

GNATURE

REORDER ALL SUPPLIES FROM N.P.A.S. P.O. BOX 1435, MANSFIELD, OH 44901

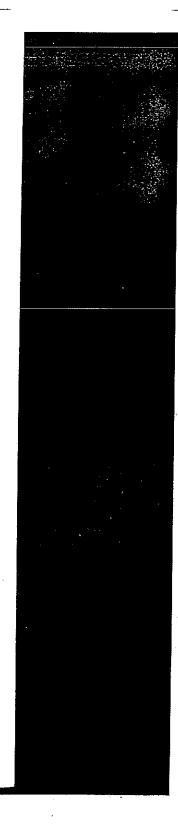

State of Missouri DEPARTMENT OF HEALTH



YPE II



Lab. 4 (R7-88)

SCOTT GARY

is hereby authorized to instruct and supervise operators, train instructors, inspect, callbrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER; ALCO-SENSOR IV/RBT IV for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986. John & Mathewson Date 03/06/09 Director of State Public Health Laboratory Number 920051 03/06/2011 Expires Director, Department of Health MO 580-0771 (7-88)

State of Missouri DEPARTMENT OF HEALTH



MO 580-0771 (7-88)

PERMIT TYPE II



Lab. 4 (R7-88)

CHRISTOPHER MCBANE

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER; ALCO-SENSOR IV

		ns 577.020 through 577.041, RSMo 1986.
Date	09/18//08	John of Mathewson
	820270	Gloster of Cialo Dubli- Usella Laboratory
Number	09/18/2010	and and
Expires		Director, Department of Health