



MISSOURI DEPARTMENT OF HEALTH
STATE PUBLIC HEALTH LABORATORY
DATAMASTER MAINTENANCE REPORT

BREATH ALCOHOL PROGRAM

JUN -3 2009

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.

DATAMASTER SN 203043	DATE OF INSPECTION 05-31-09
LOCATION OF INSTRUMENT (STREET AND CITY) Odessa MO, 101 N. 2nd St. Odessa MO 64076	TIME OF INSPECTION 11:05

CHECKLIST: Place a check (✓) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> TIME AND DATE 11:06 5/31/09	
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34 °C ± 0.2°C) 34 °C	
<input checked="" type="checkbox"/> CALIBRATION CHECK - Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE RECIRCULATION PUMP) <input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE <input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	

TEST 1 0.097	TEST 2 0.098	TEST 3 0.098
---------------------	---------------------	---------------------

<input checked="" type="checkbox"/> PERFORM R.F.I. TEST (PRINTOUT ATTACHED)							
<input checked="" type="checkbox"/> NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)							
REFUSALS 0	(0-.04) 0	(.05-.09) 0	(.10-.14) 0	(.15-.19) 0	(Over .19) 0		

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary)

Instrument operating within DHSS specs.
Repro Marketing Inc., .10%, Lot # 08001, exp. 08-11-09

INSPECTING OFFICER	
SIGNATURE Scott Ruckler	PRINT NAME Scott Ruckler
TYPE II PERMIT NUMBER/EXPIRATION DATE 820301 10/17/10	TELEPHONE NUMBER (816) 633-7575

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 08001

EXPIRATION DATE: August 11, 2009 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:


RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 08001 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed utilizing a gas chromatograph and found to contain 0.1216 gms/dl wt/vol ethyl alcohol in aqueous solution (i.e. ethanol).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of 0.100% +/-3% when heated to 34 Degrees Celsius +/-0.2 Degrees Celsius in a simulator.

The date of manufacture for this lot number is February 12, 2008.
The expiration date for this lot number is August 11, 2009 at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

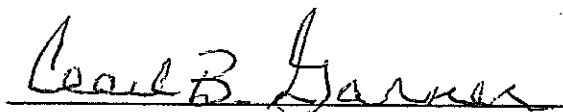

Cecil B. Garner, President
RepCo Marketing, Inc.

BAC DataMaster Evidence Ticket

STATE OF MISSOURI
ODESSA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 203043
05/31/09

TESTING OFFICER:
BUCKLER/SCOTT
OFFICER I.D.: 502
PERMIT NUMBER: 820301
EXPIRATION DATE: 10/17/10
MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---

BLANK TEST	.000	11:09
INTERNAL STANDARD	VERIFIED	11:10
EXTERNAL STANDARD	.097	11:10
BLANK TEST	.000	11:11
EXTERNAL STANDARD	.098	11:11
BLANK TEST	.000	11:12
EXTERNAL STANDARD	.098	11:12
BLANK TEST	.000	11:13

N = 3
SIM. = .1
AVG. = .0976

BAC DataMaster Evidence Ticket

STATE OF MISSOURI
ODESSA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 203043
05/31/09
11:06

--- DIAGNOSTIC CHECK ---

COMPUTER:	OKAY
PROGRAM:	OKAY
HEATERS	
SAMPLE CHAMBER:	49c
FLOW DETECTOR:	OKAY
PUMP	
HIGH SPEED:	OKAY
DETECTOR:	OKAY
FILTERS:	OKAY
QUARTZ STANDARD:	OKAY
CALIBRATION:	OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGH
IJKLMNOPQRSTUVWXYZ[\]^_`abcdefgghijklmno
pqrstuvwxyz{|}~

STATE OF MISSOURI
ODESSA POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 203043
05/31/09

ARREST TIME: 10:45
SUBJECT NAME:
RFI

DOB: 11/22/33 SEX: M

STATE/D.O.B.: MO/1234

ARRESTING OFFICER:

BUCKLER/SCOTT

OFFICER I.D.: 502

TESTING OFFICER:

BUCKLER

OFFICER I.D.: 502

PERMIT NUMBER: 820301

EXPIRATION DATE: 10/17/10

MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

BLANK TEST	.000	11:16
INTERNAL STANDARD	VERIFIED	11:16
RADIO INTERFERENCE		

Operator Signature

BAC DataMaster
Evidence Ticket

State of Missouri
DEPARTMENT OF HEALTH



P E R M I T
TYPE II



SCOTT A BUCKLER

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.

Date 10/17/08
Number 820301
Expires 10/17/2010

John J Mathewson

Director of State Public Health Laboratory

Director, Department of Health