



MISSOURI DEPARTMENT OF HEALTH
STATE PUBLIC HEALTH LABORATORY
DATAMASTER MAINTENANCE REPORT

RECEIVED Drugs Breath Alcohol Program
By Carol Day at 8:00 am, Mar 15, 2010

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.

DATAMASTER SN 204063 / INV# 127190	DATE OF INSPECTION 3-5-2010
LOCATION OF INSTRUMENT (STREET AND CITY) MARIES COUNTY SHERIFF'S DEPT. VIENNA, MO	TIME OF INSPECTION 19:37

CHECKLIST: Place a check (✓) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> TIME AND DATE 3-5-2010 @ 19:37	
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34 °C ± 0.2°C) 34.2 OK, S.N. (6-1117)	
<input checked="" type="checkbox"/> CALIBRATION CHECK - Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE RECIRCULATION PUMP) <input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE <input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	

TEST 1 0.097	TEST 2 0.096	TEST 3 0.097
---------------------	---------------------	---------------------

<input checked="" type="checkbox"/> PERFORM R.F.I. TEST (PRINTOUT ATTACHED)					
<input checked="" type="checkbox"/> NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)					
REFUSALS 0	(0-.04) 0	(.05-.09) 0	(.10-.14) 0	(.15-.19) 2	(Over .19) 0

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary)

REPLD MARKETING INC - 10% SOLUTION. LOT# 09001, EXP. DATE: 04-22-11, BOTTLE # 0137.

THIS INSTRUMENT OPERATES WITHIN THE DEPT. OF HEALTH, SENIOR SERVICES, SPECS.

INSPECTING OFFICER

SIGNATURE 	PRINT NAME B. A. VOGT, TRAINER # 987
TYPE II PERMIT NUMBER/EXPIRATION DATE 920091	TELEPHONE NUMBER (573) 368-2345

REPCO MARKETING INC.3101 JEBB STONYBROOK DRIVE
RALLIGH, NC 27604
810-870-0460

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 09001

EXPIRATION DATE: April 22, 2011 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:


RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 09001 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1214 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is April 23, 2009.
The expiration date for this lot number is April 22, 2011
at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.

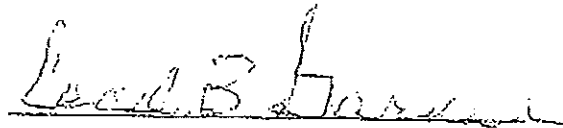

Cecil B. Garner, President
RepCo Marketing, Inc.

State of Missouri
DEPARTMENT OF HEALTH



P E R M I T
TYPE II



BROCK A VOGT

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):

DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.

Date 04/22/09
Number 920091
Expires 04/22/2011

John J Mathewson

Director of State Public Health Laboratory

Margaret T. Donnelly

Director, Department of Health

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204063
03/05/10
19:37


----- DIAGNOSTIC CHECK -----

COMPUTER: OKAY
PROGRAM (04-07-2009): OKAY
HEATERS:
SAMPLE CHAMBER: 50c
FLOW DETECTOR: OKAY
PUMP
HIGH SPEED: OKAY
DETECTOR: OKAY
FILTERS: OKAY
QUARTZ STANDARD: OKAY
CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGHIJ
KLMNOPQRSTUVWXYZ[\]^_`ab cde fgh i j k l m n o
p q r s t u v w x y z { | } ~ "

Operator Signature

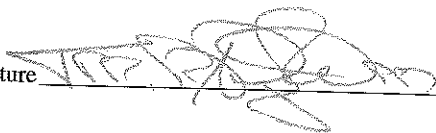

22

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204063
03/05/10


TESTING OFFICER:
VOGT/BROCK/H
OFFICER I.D.# 997
PERMIT NUMBER: 920091
EXPIRATION DATE: 04/22/11
MISCELLANEOUS DATA:

----- SUPERVISOR: NONE -----

BLANK TEST	.000	19:45
INTERNAL STANDARD	VERIFIED	19:45
EXTERNAL STANDARD	.097	19:45
BLANK TEST	.000	19:46
EXTERNAL STANDARD	.096	19:47
BLANK TEST	.000	19:47
EXTERNAL STANDARD	.097	19:48
BLANK TEST	.000	19:49

N = 3
S.D. = .1
AVG. = .0966

Operator Signature

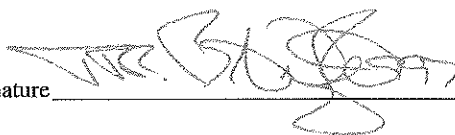

2208-02

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204063
03/05/10
19:37

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM (04-07-2009): OKAY

HEATERS

SAMPLE CHAMBER: 50c

FLOW DETECTOR: OKAY

PUMP

HIGH SPEED: OKAY

DETECTOR: OKAY

FILTERS: OKAY

QUARTZ STANDARD: OKAY

CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@ABCDEFGHIJ
KLMNOPQRSTUVWXYZ[\]^_`abcdefghijklmnopqrstuvwxyz{|}~
printuwwxyz{|}~

Operator Signature