



MISSOURI DEPARTMENT OF HEALTH
STATE PUBLIC HEALTH LABORATORY
DATAMASTER MAINTENANCE REPORT

RECEIVED
DHSS Breath Alcohol Program
By Carol Day at 4:03 pm, Mar 03, 2010

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.

DATAMASTER S.N. 204112	DATE OF INSPECTION 03-01-2010
LOCATION OF INSTRUMENT (STREET AND CITY) BROOKFIELD POLICE DEPARTMENT 116 W BROOKS, BROOKFIELD, MO	TIME OF INSPECTION 1811

CHECKLIST: Place a check (✓) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED)	
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 50 °C	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> TIME AND DATE	
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34 °C ± 0.2 °C) 34 °C	
<input checked="" type="checkbox"/> CALIBRATION CHECK - Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE RECIRCULATION PUMP) <input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE <input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	

TEST 1 .099	TEST 2 .099	TEST 3 .100
-------------	-------------	-------------

<input checked="" type="checkbox"/> PERFORM R.F.I. TEST (PRINTOUT ATTACHED)						
<input checked="" type="checkbox"/> NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)						
REFUSALS 0	(0-.04) 0	(.05-.09) 0	(.10-.14) 0	(.15-.19) 1	(Over .19) 0	

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary)

MEETS DEPARTMENT OF HEALTH STANDARDS

REPCO MARKETING, INC
VAPOR SOLUTION → .100

LOT # 09002
EXP: 08-31-2011

INSPECTING OFFICER	
SIGNATURE Reid T Fessler	PRINT NAME REID T FESSLER
TYPE II PERMIT NUMBER/EXPIRATION DATE 200032 02-05-2012	TELEPHONE NUMBER 660-385-2132

CERTIFICATE OF ANALYSIS

MANUFACTURER AND SUPPLIER: RepCo Marketing, Inc.

LOT NUMBER: 09002

EXPIRATION DATE: August 31, 2011 at 11:59 p.m.

RepCo Marketing, Inc. certifies the following:

RepCo Marketing, Inc. manufactured, tested and supplied Lot Number 09002 of Alcohol Certified Solution for simulators. Random samples of said lot number were analyzed by an independent laboratory utilizing a gas chromatograph and found to contain .1206 gms/dl +/- .003 gms/dl wt/vol ethanol (95% Confidence).

The alcohol and distilled water used in the solution were found to be free of any interfering substance.

This solution will produce a vapor alcohol value of .100 +/- 3% gms/210L Breath when heated to 34 Degrees Celsius +/- 0.2 Degrees Celsius in a simulator (95% Confidence).

The date of manufacture for this lot number is September 1, 2009.
The expiration date for this lot number is August 31, 2011
at 11:59 p.m.

This document is a true representation of the original Certificate of Analysis.



Cecil B. Garner, President
RepCo Marketing, Inc.

State of Missouri
DEPARTMENT OF HEALTH



P E R M I T
TYPE II




REID T FESSLER

is hereby authorized to instruct and supervise operators, train instructors, inspect, calibrate, perform field repairs, and operate the following breath analyzer(s):


DATAMASTER

for the determination of the alcoholic content of blood from a sample of expired (alveolar) air. Issued under the provisions of sections 577.020 through 577.041, RSMo 1986.

Date 02/05/2010
Number 200032
Expires 02/05/2012

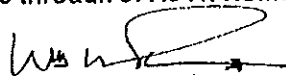 Interim Director

Director of State Public Health Laboratory

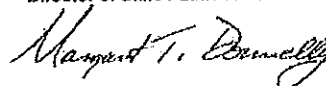

Director, Department of Health

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
03/01/10

ARREST TIME: 08:00

SUBJECT NAME:

TEST

DOB: 11/11/11 SEX: M

STATE/D.L.: MO/TEST

ARRESTING OFFICER:

FESSLER/R/T

OFFICER I.D.: 992

TESTING OFFICER:

FESSLER/R/T

OFFICER I.D.: 992

PERMIT NUMBER: 200002

EXPIRATION DATE: 02/05/12

MISCELLANEOUS DATA:

--- BREATH ANALYSIS ---

RADIO INTERFERENCE

Operator Signature

R.T. Fessler

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
03/01/10
18:11

--- DIAGNOSTIC CHECK ---

COMPUTER: OKAY

PROGRAM (04-07-2003): OKAY

HEATERS

SAMPLE CHAMBER: 50c

FLOW DETECTOR: OKAY

PUMP

HIGH SPEED: OKAY

DETECTOR: OKAY

FILTERS: OKAY

QUARTZ STANDARD: OKAY

CALIBRATION: OKAY

PRINTER TEST

!"#\$%&'()*+,-./0123456789:;<=>?@A8CDEFG
HIJKLMNOPQRSTUVWXYZ[\]^_`abcedefghijklmno
pqrs tuvxyz{|}~"

Operator Signature

R.T. Fessler

Face This Side Down - This Edge In First

BAC DataMaster Evidence Ticket

MISSOURI STATE HIGHWAY PATROL
BAC DATAMASTER SERIAL NUMBER 204112
03/01/10

TESTING OFFICER:
FESSLER/R/T
OFFICER I.D.: 992
PERMIT NUMBER: 200032
EXPIRATION DATE: 02/05/12
MISCELLANEOUS DATA:
.10 TEST

--- SUPERVISOR NODE ---

BLANK TEST	.000	18:15
INTERNAL STANDARD	VERIFIED	18:15
EXTERNAL STANDARD	.099	18:16
BLANK TEST	.000	18:16
EXTERNAL STANDARD	.099	18:17
BLANK TEST	.000	18:17
EXTERNAL STANDARD	.100	18:18
BLANK TEST	.000	18:18

N = 3

STDEV. = .1

AVG. = .0999

Operator Signature

