



MISSOURI DEPARTMENT OF HEALTH  
STATE PUBLIC HEALTH LABORATORY  
**DATAMASTER MAINTENANCE REPORT**

RECEIVED  
Missouri Breath Alcohol Program  
By Carol Day at 8:15 am, Mar 04, 2010

Complete this report in duplicate at the time of the regular monthly preventive maintenance check, and whenever instrument is repaired. Send copy to Department of Health; retain original in department file.

DATAMASTER SN #940140	DATE OF INSPECTION 01-15-10
LOCATION OF INSTRUMENT (STREET AND CITY) 1414 EVERGREEN PARK - WELSTON	TIME OF INSPECTION 2330

**CHECKLIST:** Place a check (✓) to the left of each item if found to be satisfactory or if operating within established limits. (Write in observed values where determined.) Unchecked items must be corrected before using instrument.

<input checked="" type="checkbox"/> DIAGNOSTIC CHECK (PRINTOUT ATTACHED) PASS	
<input checked="" type="checkbox"/> COMPUTER	<input checked="" type="checkbox"/> DETECTOR
<input checked="" type="checkbox"/> PROGRAM	<input checked="" type="checkbox"/> FILTERS
<input checked="" type="checkbox"/> HEATERS SAMPLE CHAMBER 49 °C PASS	<input checked="" type="checkbox"/> QUARTZ STANDARD
<input checked="" type="checkbox"/> FLOW DETECTOR	<input checked="" type="checkbox"/> CALIBRATION
<input checked="" type="checkbox"/> PUMP HIGH SPEED	<input checked="" type="checkbox"/> PRINTER
<input checked="" type="checkbox"/> INDICATOR LIGHTS	
<input checked="" type="checkbox"/> TIME AND DATE	
<input checked="" type="checkbox"/> SIMULATOR TEMPERATURE (34 °C ± 0.2°C) 34°C PASS	
<input checked="" type="checkbox"/> CALIBRATION CHECK - PASS Run three tests using a standard solution. All three tests must be within ± 5% of the standard value and must have a spread of .005 or less. Check the box corresponding to the standard solution being used. (PRINTOUT ATTACHED) (USE RECIRCULATION PUMP) <input checked="" type="checkbox"/> 0.100% STANDARD - MUST READ BETWEEN 0.095% and 0.105% INCLUSIVE <input type="checkbox"/> 0.040% STANDARD - MUST READ BETWEEN 0.038% and 0.042% INCLUSIVE (ONLY ONE STANDARD IS TO BE USED PER MAINTENANCE REPORT)	

TEST 1 ✓ .095	TEST 2 ✓ .096	TEST 3 ✓ .095
---------------	---------------	---------------

☒ PERFORM R.F.I. TEST (PRINTOUT ATTACHED) PASS

☒ NUMBER OF REFUSALS, SINCE LAST MAINTENANCE REPORT, AND NUMBER OF BREATH TESTS IN EACH RANGE AS FOLLOWS: (DO NOT INCLUDE SIMULATOR TESTS)

REFUSALS 1	(0-.04) 0	(.05-.09) 0	(.10-.14) 1	(.15-.19) 0	(Over .19) 1
------------	-----------	-------------	-------------	-------------	--------------

List any new parts and describe any alteration or modification that was made to restore the instrument to operate satisfactorily and within established limits (use other side if necessary)

INSTRUMENT WITHIN D.O.H. GUIDELINES  
RESET PRINTER REBORN

GUTH LABS lot # 09270 RES. 10 SALVAGE EXP - 09-23-10

<b>INSPECTING OFFICER</b>	
SIGNATURE P. J. BRUGMANN 662	PRINT NAME P. J. BRUGMANN 662
TYPE II PERMIT NUMBER/EXPIRATION DATE #920015 02-06-11	TELEPHONE NUMBER (314) 553-8010



**GUTH LABORATORIES, INC.**

590 NORTH 67th STREET • HARRISBURG, PA 17111-4511 • TELEPHONE: 717-564-5470

## CERTIFICATE OF ANALYSIS

Certified Alcohol Reference Solution for Simulator

Random Samples of Lot Number **09270** of Alcohol Reference Solution for Simulator were analyzed by gas chromatography and found to contain **0.1207** percent (w/vol) ethyl alcohol. The expiration date for this lot number is **September 23, 2010** at 11:59 PM.

When used in a calibrated Simulator, operating at  $34^{\circ}\text{C} \pm .2^{\circ}\text{C}$ , this solution will give a breath alcohol analysis instrument reading of **0.10** percent BAC.

The alcohol and water used in this solution were free of test interfering substances.

Ted L. Pauley, President  
GUTH LABORATORIES, INC.

Face This Side Down - This Edge In First

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
WELLSTON POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 940140  
01/15/10

ARREST TIME: 23:08

SUBJECT NAME:

TEST

DOB: 01/01/01 SEX: M

STATE/D.L.: MO/TEST

ARRESTING OFFICER:

TEST

OFFICER I.D.: TEST

TESTING OFFICER:

BRUGHMAN

OFFICER I.D.: 662

PERMIT NUMBER: 920015

EXPIRATION DATE: 02/06/11


MISCELLANEOUS DATA:

RFI TEST

--- BREATH ANALYSIS ---

BLANK TEST	.000	23:06
INTERNAL STANDARD	VERIFIED	23:06
RADIO INTERFERENCE		

Operator Signature

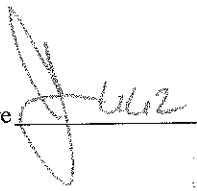

2208-02

Face This Side Down - This Edge In First

# BAC DataMaster Evidence Ticket

STATE OF MISSOURI  
WELLSTON POLICE DEPARTMENT

BAC DATAMASTER SERIAL NUMBER 940140  
01/15/10

TESTING OFFICER:

BRUGHMAN

OFFICER I.D.: 662

PERMIT NUMBER: 920015

EXPIRATION DATE: 02/06/11

MISCELLANEOUS DATA:

--- SUPERVISOR MODE ---


BLANK TEST	.000	23:51
INTERNAL STANDARD	VERIFIED	23:51
EXTERNAL STANDARD	.000	23:51
BLANK TEST	.000	23:52
EXTERNAL STANDARD	.000	23:53
BLANK TEST	.000	23:53
EXTERNAL STANDARD	.000	23:54
BLANK TEST	.000	23:55

N = 2

SD = .1

AVG = .0950

Operator Signature

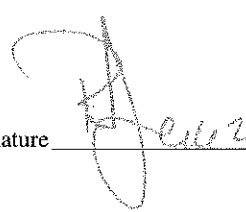

2208-02

# BAC DataMaster

## Evidence Ticket

BRF INSTRUCTOR SERIAL NUMBER 340140  
01/15/10  
23:30

... .. DIAGNOSTIC CHECK ... ..

**CONCLUSION**

```

1"=32K"04+-, 3123456789101112131415161718192021222324252627282930313233343536373839404142434445464748495051525354555657585960616263646566676869707172737475767778798081828384858687888990919293949596979899100101102103104105106107108109110111112113114115116117118119120121122123124125126127128129130131132133134135136137138139140141142143144145146147148149150151152153154155156157158159160161162163164165166167168169170171172173174175176177178179180181182183184185186187188189190191192193194195196197198199200201202203204205206207208209210211212213214215216217218219220221222223224225226227228229230231232233234235236237238239240241242243244245246247248249250251252253254255256257258259260261262263264265266267268269270271272273274275276277278279280281282283284285286287288289290291292293294295296297298299300301302303304305306307308309310311312313314315316317318319320321322323324325326327328329330331332333334335336337338339340341342343344345346347348349350351352353354355356357358359360361362363364365366367368369370371372373374375376377378379380381382383384385386387388389390391392393394395396397398399400401402403404405406407408409410411412413414415416417418419420421422423424425426427428429430431432433434435436437438439440441442443444445446447448449450451452453454455456457458459460461462463464465466467468469470471472473474475476477478479480481482483484485486487488489490491492493494495496497498499500501502503504505506507508509510511512513514515516517518519520521522523524525526527528529530531532533534535536537538539540541542543544545546547548549550551552553554555556557558559560561562563564565566567568569570571572573574575576577578579580581582583584585586587588589590591592593594595596597598599600601602603604605606607608609610611612613614615616617618619620621622623624625626627628629630631632633634635636637638639640641642643644645646647648649650651652653654655656657658659660661662663664665666667668669670671672673674675676677678679680681682683684685686687688689690691692693694695696697698699700701702703704705706707708709710711712713714715716717718719720721722723724725726727728729730731732733734735736737738739740741742743744745746747748749750751752753754755756757758759760761762763764765766767768769770771772773774775776777778779780781782783784785786787788789790791792793794795796797798799800801802803804805806807808809810811812813814815816817818819820821822823824825826827828829830831832833834835836837838839840841842843844845846847848849850851852853854855856857858859860861862863864865866867868869870871872873874875876877878879880881882883884885886887888889890891892893894895896897898899900901902903904905906907908909910911912913914915916917918919920921922923924925926927928929930931932933934935936937938939940941942943944945946947948949950951952953954955956957958959960961962963964965966967968969970971972973974975976977978979980981982983984985986987988989990991992993994995996997998999100010011002100310041005100610071008100910101011101210131014101510161017101810191020102110221023102410251026102710281029103010311032103310341035103610371038103910401041104210431044104510461047104810491050105110521053105410551056105710581059106010611062106310641065106610671068106910701071107210731074107510761077107810791080108110821083108410851086108710881089109010911092109310941095109610971098109911001101110211031104110511061107110811091110111111121113111411151116111711181119112011211122112311241125112611271128112911301131113211331134113511361137113811391140114111421143114411451146114711481149115011511152115311541155115611571158115911601161116211631164116511661167116811691170117111721173117411751176117711781179118011811182118311841185118611871188118911901191119211931194119511961197119811991200120112021203120412051206120712081209121012111212121312141215121612171218121912201221122212231224122512261227122812291230123112321233123412351236123712381239124012411242124312441245124612471248124912501251125212531254125512561257125812591260126112621263126412651266126712681269127012711272127312741275127612771278127912801281128212831284128512861287128812891290129112921293129412951296129
```

Operator Signature